
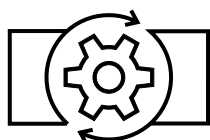


ChimerisMonitor IVD


ChimerisMonitor RUO

The intuitive software solutions for automated chimerism analysis

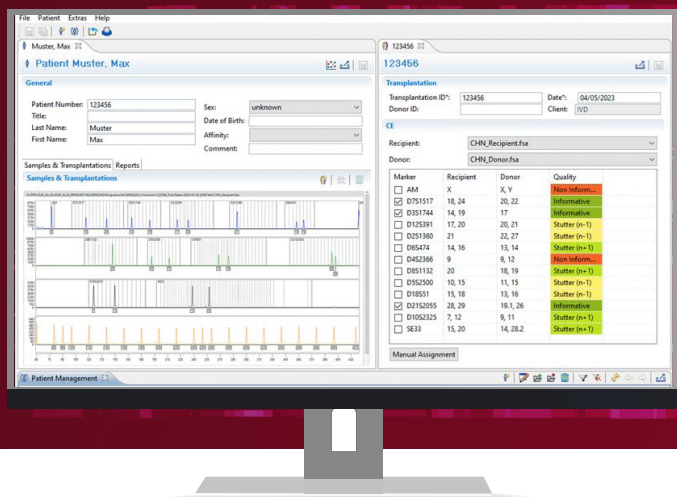


Easy

Simplifies and speeds up data analysis for efficient BIOTYPE chimerism kit evaluation

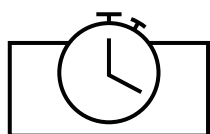

Fast

Automates critical and time-consuming steps from allelic call to chimerism value calculation

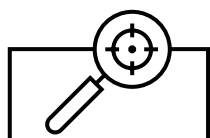

Accurate

Extensive control concept for run and sample validation



Comprehensive

Builds a patient database encompassing preTx genotypes and postTx chimerism kinetics

ChimerisMonitor product details

PRODUCT DESCRIPTION

ChimerisMonitor is the professional software for automated chimerism analysis combining speed and accuracy with an intuitive and customized interface. The software mediates analysis of chimerism samples from initial genotyping to highly sensitive quantification of residual donor or recipient cells. It automatically identifies donor or recipient signals and provides a ranking of identified loci in regard to their informativity. It calculates percent-chimerism from various specimens and depicts the standard deviation of each analysis. Obtained results are summarized in a comprehensive report conferring clarity and providing overview in longitudinal studies.

FLEXIBLE SOLUTION FOR YOUR NEEDS

	 ChimerisMonitor RUO	 ChimerisMonitor CE-IVD	 ChimerisMonitor Bundle
Analysis of Mentype® Chimera® and Mentype® DIPscreen	✓	✓	✓
Analysis of Mentype® DIPquants	✓	✗	✓
IVD-certified workflow for secure, robust & reliable data evaluation	✗	✓	✓
Extensive control concept for run and sample validation	✗	✓	✓
Visualization of EPGs and Size calling regressions incl. attachment (controls) to every sample	✓	✓	✓
Patient management at one glance and summary of all monitoring data in patient-specific display including report (pdf, csv)	✓	✓	✓
Synchronized database including visualization of data status (IVD/RUO)	✗	✗	✓
Reanalysis of critical data sets & reassignment of EPGs	✓	✗	✓

ORDER INFORMATION

Product	Size	Status
ChimerisMonitor IVD	46-14800-0000	CE-IVD
ChimerisMonitor RUO	46-14801-0000	RUO

Direct your orders via email to sales@biotype.de



Scan the code to learn more about ChimerisMonitor software